

# DeVilbiss IntelliPAP Series: Single Battery Connection


DeVilbiss IntelliPAP series devices can be connected directly to the Freedom™ CPAP Battery using the DC power pigtail cord and appropriate connector tip included in the Freedom™ CPAP Battery Kit. Please note IntelliPAP 2 series devices will require the use of the 150W sine wave power inverter in conjunction with the Freedom™ CPAP Battery Kit to complete the connection between the PAP device and the battery.

## REQUIRED COMPONENTS

- Freedom™ CPAP Battery Kit

## CONNECTION INSTRUCTIONS

- Plug the DC power pigtail cord into the DC OUT jack of the battery.
- Insert the connector tip with blue top labeled “J” into the DC power pigtail cord.
- Plug the DC power pigtail cord with “J” connector tip attached into the DC input port of the PAP device.
- Turn the PAP device on and it should begin functioning.
- The digital display screen will show the charge level of the battery. You may press the charge level indicator button to illuminate this screen. It will stay illuminated for a few seconds at a time.
- When not in use disconnect all devices from the battery thus allowing the battery to completely shut down.

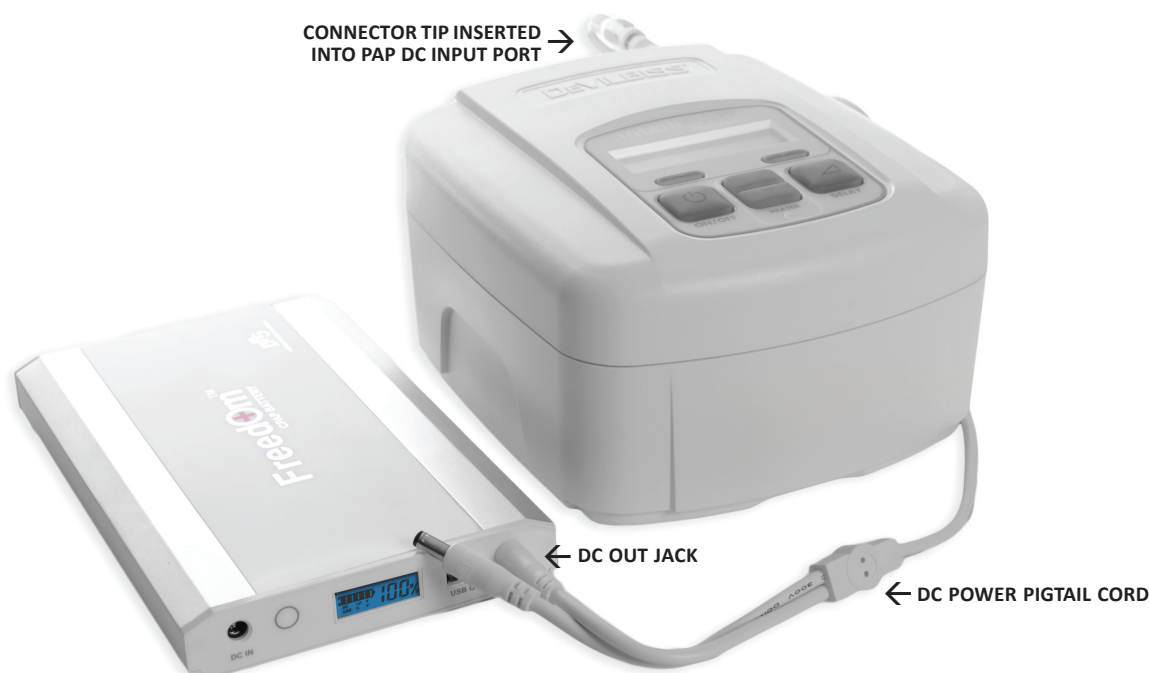

## OPTIONAL FREEDOM™ CPAP BATTERY CHARGING ACCESSORIES

- Auto Charge DC Cable {BPS DC.CABLE}: Used to charge the battery in a vehicle (vehicle must be running).
- 150W Sine Wave Power Inverter {BPS 150 INVERTER}: Used to charge the battery in a vehicle (vehicle does not need to be running) and allows the battery to power/charge other small electronic devices.
- SunPower 40W Solar Panel {BPS SP40}: Used to charge the battery using the power of the sun.

# DeVilbiss IntelliPAP Series: Dual Battery Connection


DeVilbiss IntelliPAP series devices can be connected directly to the Freedom™ CPAP Battery using the DC power pigtail cord and appropriate connector tip included in the Freedom™ CPAP Battery Kit. Please note IntelliPAP 2 series devices will require the use of the 150W sine wave power inverter in conjunction with the Freedom™ CPAP Battery Kit to complete the connection between the PAP device and the battery.

## REQUIRED COMPONENTS

- Freedom™ CPAP Battery Kit
- Freedom™ CPAP Battery (additional battery only)

## CONNECTION INSTRUCTIONS

- Plug the DC power pigtail cord into the DC OUT jacks of both batteries.
- Insert the connector tip with blue top labeled “J” into the DC power pigtail cord.
- Plug the DC power pigtail cord with “J” connector tip attached into the DC input port of the PAP device.
- Turn the PAP device on and it should begin functioning.
- The digital display screen will show the charge level of the battery. You may press the charge level indicator button to illuminate this screen. It will stay illuminated for a few seconds at a time.
- When not in use disconnect all devices from the battery thus allowing the battery to completely shut down.

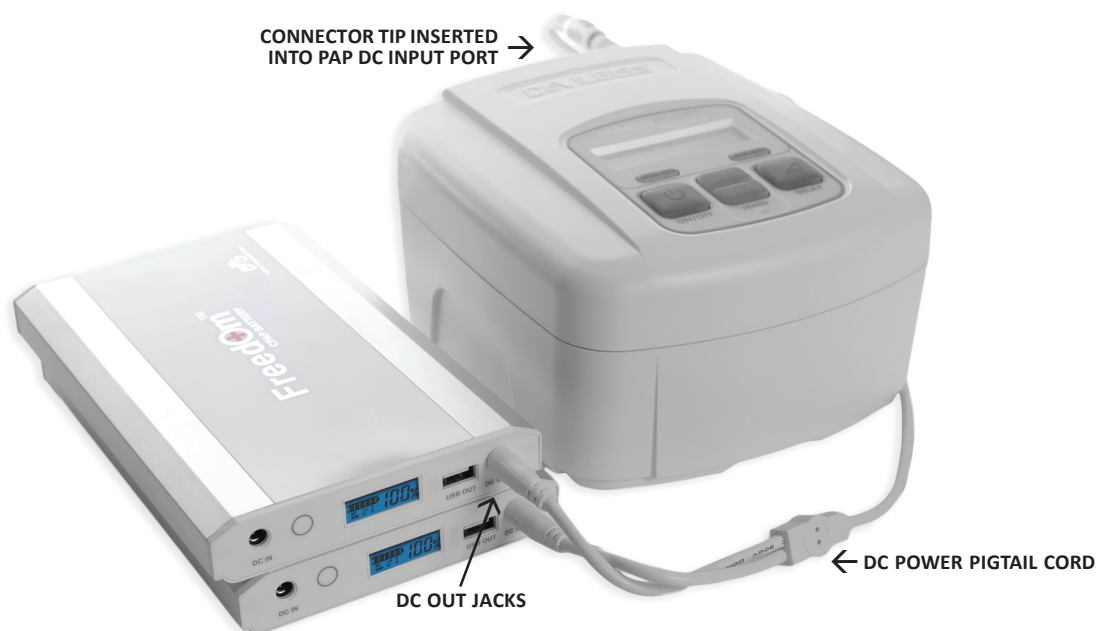

## OPTIONAL FREEDOM™ CPAP BATTERY CHARGING ACCESSORIES

- Auto Charge DC Cable {BPS DC.CABLE}: Used to charge the battery in a vehicle (vehicle must be running).
- 150W Sine Wave Power Inverter {BPS 150 INVERTER}: Used to charge the battery in a vehicle (vehicle does not need to be running) and allows the battery to power/charge other small electronic devices.
- SunPower 40W Solar Panel {BPS SP40}: Used to charge the battery using the power of the sun.