

# Respironics M & System One Series: Single Battery Connection

Respironics System One and M series devices can be connected directly to the Freedom™ CPAP Battery using the DC power pigtail cord and the appropriate connector tip included in the Freedom™ CPAP Battery Kit.


## REQUIRED COMPONENTS

- Freedom™ CPAP Battery Kit

## CONNECTION INSTRUCTIONS

- Plug the DC power pigtail cord into the DC OUT jack of the battery.
- Insert the appropriate connector tip into the DC power pigtail cord.
  - System One 60 Series: Tip with black top labeled “W”
  - System One 50 Series & M Series: Tip with yellow top labeled “N”
- Plug the DC power pigtail cord with appropriate connector tip attached into the DC input port of the PAP device.
- Turn the PAP device on and it should begin functioning.
- The digital display screen will show the charge level of the battery. You may press the charge level indicator button to illuminate this screen. It will stay illuminated for a few seconds at a time.
- When not in use disconnect all devices from the battery thus allowing the battery to completely shut down.

Image shows connection to Respironics System One 60 series PAP device.

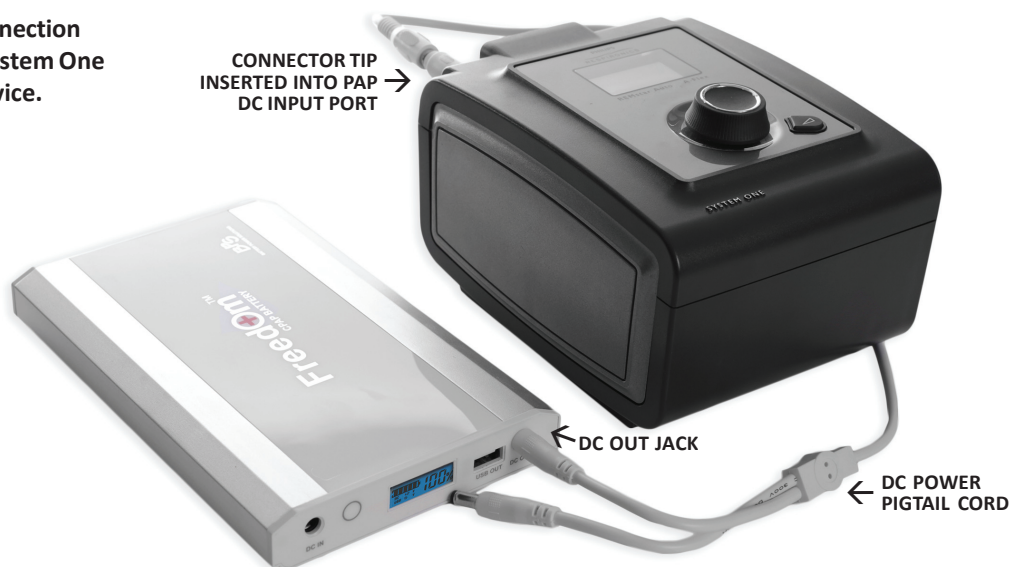

## OPTIONAL FREEDOM™ CPAP BATTERY CHARGING ACCESSORIES

- Auto Charge DC Cable {BPS DC.CABLE}: Used to charge the battery in a vehicle (vehicle must be running).
- 150W Sine Wave Power Inverter {BPS 150 INVERTER}: Used to charge the battery in a vehicle (vehicle does not need to be running) and allows the battery to power/charge other small electronic devices.
- SunPower 40W Solar Panel {BPS SP40}: Used to charge the battery using the power of the sun.

# Respironics M & System One Series: Dual Battery Connection

Respironics System One and M series devices can be connected directly to the Freedom™ CPAP Battery using the DC power pigtail cord and the appropriate connector tip included in the Freedom™ CPAP Battery Kit.


## REQUIRED COMPONENTS

- Freedom™ CPAP Battery Kit
- Freedom™ CPAP Battery (additional battery only)

## CONNECTION INSTRUCTIONS

- Plug the DC power pigtail cord into the DC OUT jacks of both batteries.
- Insert the appropriate connector tip into the DC power pigtail cord.
  - System One 60 Series: Tip with black top labeled “W”
  - System One 50 & M Series: Tip with yellow top labeled “N”
- Plug the DC power pigtail cord with appropriate connector tip attached into the DC input port of the PAP device.
- Turn the PAP device on and it should begin functioning.
- The digital display screen will show the charge level of the battery. You may press the charge level indicator button to illuminate this screen. It will stay illuminated for a few seconds at a time.
- When not in use disconnect all devices from the battery thus allowing the battery to completely shut down.

Image shows connection to Respironics System One 60 series PAP device.

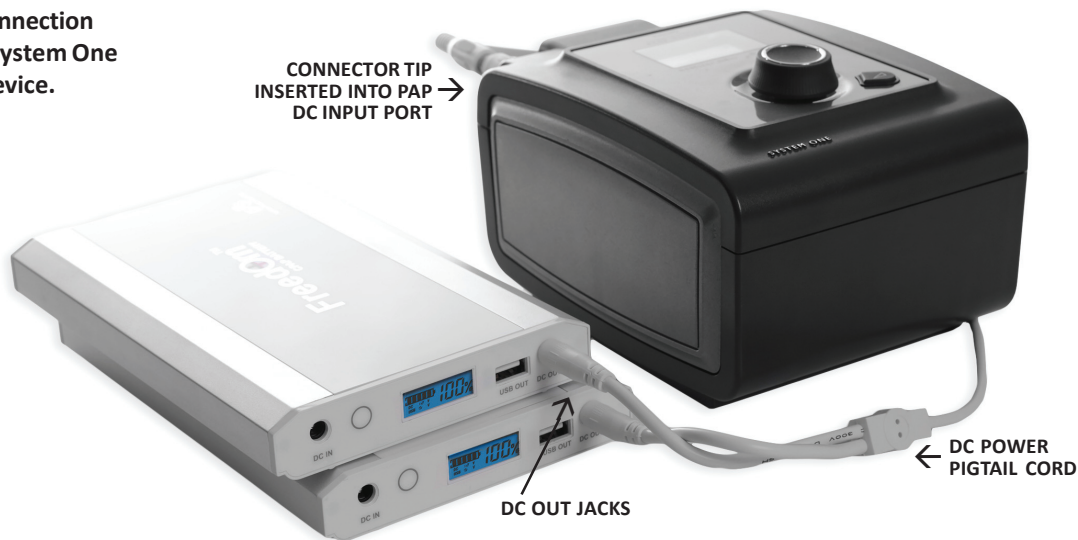

## OPTIONAL FREEDOM™ CPAP BATTERY CHARGING ACCESSORIES

- Auto Charge DC Cable {BPS DC.CABLE}: Used to charge the battery in a vehicle (vehicle must be running).
- 150W Sine Wave Power Inverter {BPS 150 INVERTER}: Used to charge the battery in a vehicle (vehicle does not need to be running) and allows the battery to power/charge other small electronic devices.
- SunPower 40W Solar Panel {BPS SP40}: Used to charge the battery using the power of the sun.